


F35 R²



F35 R²

Compact folding wheelchair with great portability and increased range

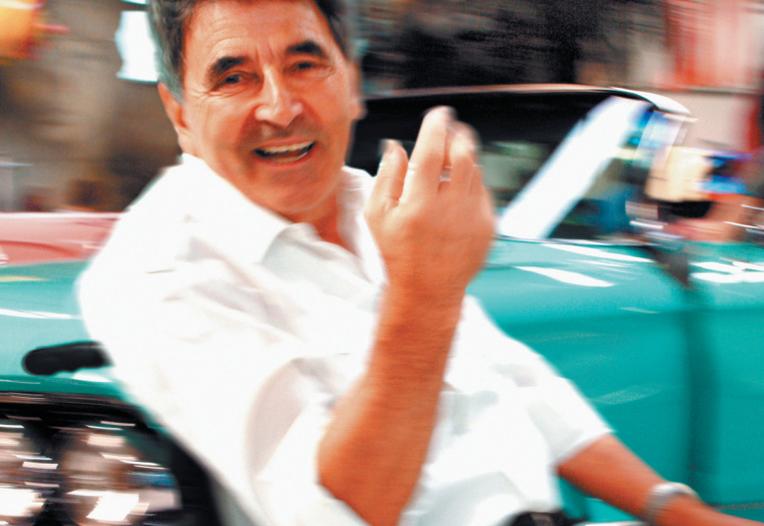

Stylish, comfortable and folding power wheelchair

The F35R² has multiple features to provide maximum comfort and accessibility in your life.

50ah batteries combined with powerful motors offer enhanced range and performance.

Easily Transportable

After removing the battery boxes and detaching the backrest by using a simple lever design, the chair can be folded to assist in transportation and storage.

The anti-tip tubes are also detachable and make the chair compact to stow away in the boot of the car.



Compact folding and easy to transport after removing the battery boxes and the backrest.

Tailored power chair

The F35R² is tailored to adapt to the needs of each user. The armrests provide a width adjustment of 4,5 cm each side.

This allows the standard seat width to have 9 cm of adjustment to a possible seat width of 52 cm. Height adjustment is achieved using a simple lever mechanism.




Width adjustable armrests, allow adjustment up to 52 cm



Height adjustable armrests, ensure maximum comfort



The footrests can swing in and out

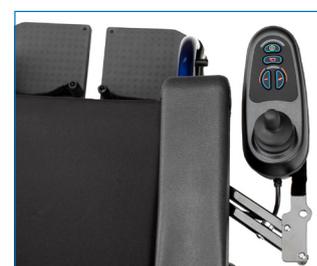

Parallel swing-away control box.



Ergonomic designed armrest

Go everywhere...

The footrests are detachable and can swing either inwards or outwards. Thanks to this feature, the length of the chair can be altered without removing the footrest prior to entering small spaces such as lifts.

Seat cushion is included to give more comfort.



Reclining backrest

Available options

- Kerb climber
- Attendant control only
- Dual control (user and attendant)
- Chin control
- Elevating footrest
- Reclining backrest
- Solid rear and front wheels

Colour range



Blue



Gold




Summer red

F35 R²

Compact folding wheelchair with great portability and increased range

Technical specifications

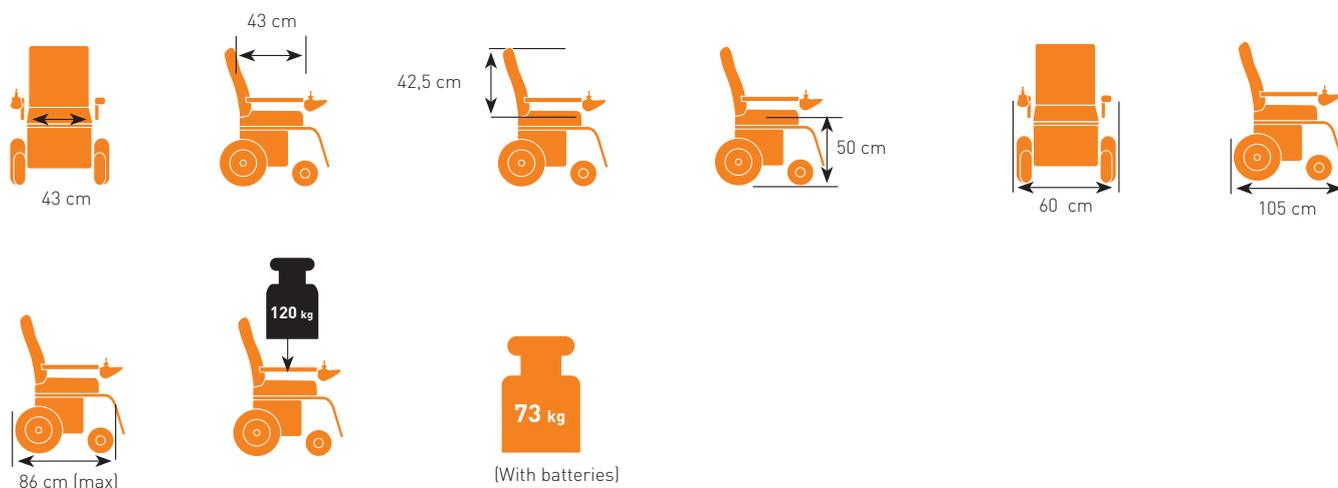

Seat width: 43 cm (can be increased up to 52cm with the armrest)	Overall width folded:	33,5cm
Seat depth: 43 cm	Turning circle:	60 cm
Backrest height: 42,5 cm	Max Slope:	10° (18%)
Height seat to floor: 50 cm	Batteries:	50 Amp
Overall width: 60 cm	Charger:	8 Ah
Overall length: 105 cm	Frame colour:	Blue, gold and summer red
Overall length without footrest: 86 cm (max)	Max. user weight:	120 kg

Tolerance +/- 1 cm

For further information on the full specification, options and accessories, please refer to the order form.
All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.
If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.com or alternatively is available on request in large text



Sunrise Medical S.L
Polígono Bakiola, 41
48498 Arrankudiaga
Vizcaya, España

Tel: 902 14 24 34
Fax: 94 648 15 75

www.SunriseMedical.es



Rev.2- July 2014